

Ammonia Vacu-vials® Kit

K-1413: 0 - 3.00 ppm N (Prog. # 12)

0 - 60.0 ppm N (Prog. # 13)

Instrument Set-up

For CHEMetrics photometers, follow the **Setup and Measurement Procedures** in the operator's manual. For spectrophotometers, set the wavelength to 610 nm. A sealed ZERO ampoule is supplied in this kit for zeroing when the sample is colorless and not turbid. For improved accuracy with colored or turbid samples, Sample Zeroing Accessory Pack, Cat. # A-0025 is recommended. Fill the A-0025 test tube with the sample (diluted sample for the 0-60 ppm range) and use this in place of the supplied ZERO ampoule to zero the instrument.

0 - 3 ppm Test Procedure

1. Fill the sample cup to the 20 mL mark with the sample to be tested. Fig. 1
2. Add 4 drops of A-1404 Stabilizer Solution. Fig. 2
3. Add 4 drops of A-1405 Catalyzer Solution (green). Fig. 2
4. Add 4 drops of A-1406 Activator Solution (blue). Fig. 2
5. **Immediately** place the Vacu-vial ampoule, tip first, into the sample cup. Stir briefly to mix the contents of the cup, then snap the tip. The ampoule will fill leaving a bubble for mixing. Fig. 3
6. To mix the ampoule, invert it several times, allowing the bubble to travel from end to end.
7. Dry the ampoule. Obtain a test result **5 minutes** after snapping the tip.

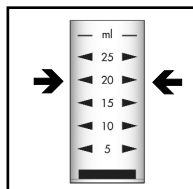


Figure 1

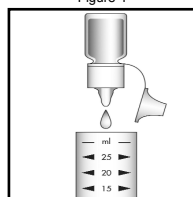


Figure 2

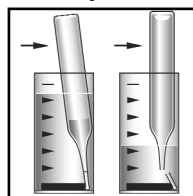


Figure 3

8. Insert the Vacu-vial ampoule into the photometer, flat end first, and obtain a reading in ppm (mg/Liter) ammonia-nitrogen (NH₃-N).

NOTE: If using a spectrophotometer that is not pre-calibrated for CHEMetrics products, then use the appropriate equation below or the **Concentration Calculator** found under the Support tab at www.chemetrics.com.

0 - 3 ppm: ppm = 2.22 (abs) - 0.07

0 - 60 ppm: ppm = 44.4 (abs) - 1.4

0 - 60 ppm Test Procedure

1. Using the syringe provided, obtain **1.0 mL** of the sample to be tested and dispense it into the empty sample cup.
2. Dilute the contents of the sample cup to the **20 mL** mark with distilled water.
3. Perform the 0 - 3 ppm Test Procedure, beginning with Step 2.

Test Method

The Ammonia Vacu-vials®¹ test kit employs the Hydroxybenzyl alcohol (HBA) chemistry.² Free ammonia reacts with hypochlorite to form monochloramine. Monochloramine reacts with HBA, in the presence of sodium nitro-ferricyanide, to form a green colored complex. This test method measures the sum of free ammonia and monochloramine. High levels of ammonia can produce false low or off color test results. Dilute the sample if the ammonia concentration is suspected to significantly exceed the test range.

1. Vacu-vials is a registered trademark of CHEMetrics, Inc. U.S. Patent No. 3,634,038
2. Krom, Michael D., Spectrophotometric Determination of Ammonia: A study of a Modified Berthelot Reduction Using Salicylate and Dichloroisocyanurate, *The Analyst*, V105 pp. 305-316, 1980.

Safety Information

Read SDS (available at www.chemetrics.com) before performing this test procedure. Wear safety glasses and protective gloves.

Visit www.chemetrics.com to view product demonstration videos.
Always follow the test procedure above to perform a test.



Simplicity in Water Analysis

www.chemetrics.com

4295 Catlett Road, Midland, VA 22728 U.S.A.

Phone: (800) 356-3072; Fax: (540) 788-4856

E-Mail: orders@chemetrics.com

Sept. 19, Rev. 3