



## Quick Operation Guide


Bench-top Smart Vacuum Aspiration System

### Smart Vac

#### Packing List

Main unit x 1, Bottle for liquid collection x 1, Tube for vacuum x 1, Tube for absorbing liquid x 1, Air filter x 2, Tip detrusion adapter x 1, 1ml Tip x 1, Quick guide x 1, Power Adapter x 1.

#### Introduction

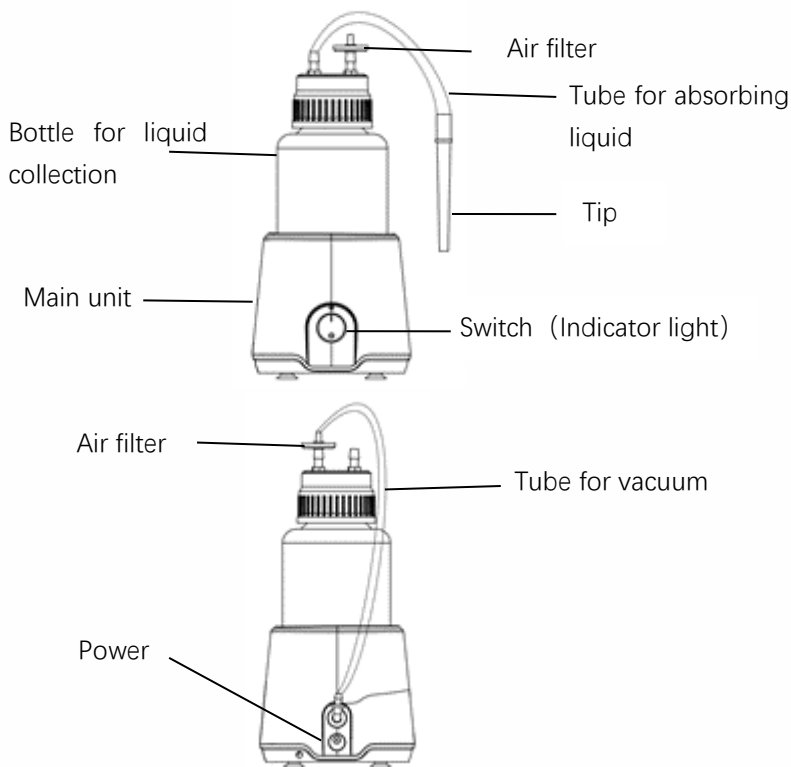



## Operation

- 1) Installation as shown above. Make sure all parts put into position.
- 2) Insert the 1ml tip into the Tip detrusion adapter.
- 3) Turn on the power switch, the indicator light is on, and start vacuuming. Insert the tip into the liquid and aspire the liquid into the bottle.
- 4) Turn off the power switch.
- 5) Take off the bottle, open the lid, and pour out the liquid.

## Notice

1. Before using for the first time, please charge with the supplied power adapter for 10 hours.
2. If the instrument is not used for a long time, please charge the battery every 2 months to avoid battery damage.
3. Please use the original power adapter for charging or plugging in. Inappropriate voltage/current will damage the internal motor.
4. After replacing the bottle, please make sure that the cap is tightened to avoid insufficient liquid absorption capacity due to air leakage.
5. It is forbidden to suck liquid into the vacuum pump of the main unit to avoid the motor damage.
6. Regularly replace the filter to keep airflow smoothly, avoid reducing the service life of the vacuum pump due to blocked pipelines.
7. When there is too much liquid in the bottle, please pour out the waste liquid in time to prevent the liquid from being sucked into the vacuum pump.

