



Forceps are metallic instruments with two branches that are usually used to catch, hold, attract or compress tissues or other surgical instruments.

Nahita dissecting forceps are made of stainless steel, and specially aimed to dissection of small animals in education

▶ Curved London-College dissecting forceps



1 | With curved tip to lift portions of tissues more freely, and striated to facilitate tissue holding.

Code	Length
FMG001	115 mm
FMG002	150 mm

▶ Laboratory forceps



1 | Curved tips, made of stainless steel

Code	Length
FMG003	110 mm

▶ Mosquito forceps, curved



1 | Made of stainless steel

Code	Length
FMG004	120 mm

▶ Dissecting forceps without teeth



1 | With straight and striated tip to facilitate tissue holding.

Code	Length
FMJ001	125 mm
FMJ002	140 mm
FMJ003	160 mm
FMJ004	200 mm
FMJ005	250 mm



▶ Dissecting forceps, thin tip



1 | Ideal for small vessels and structures.

Code	Length
FMJ006	125 mm, striated tip
FMJ007	140 mm, striated tip

▶ Dissecting forceps with teeth



1 | Instead of a striated tip, they present teeth on one side that fits into the teeth of the opposite side of the tip.
2 | These forceps are usually used to hold gauzes but they can also be used to firmly hold hard tissues and skin.

Code	Teeth	Length
FMJ008	1:2	125 mm
FMJ009	1:2	140 mm
FMJ010	1:2	160 mm
FMJ011	1:2	200 mm
FMJ012	2:3	130 mm
FMJ013	2:3	145 mm
FMJ014	1:2	140 mm*

*thin tip

▶ Adson forceps



Code	Length
FML001	120 mm

▶ Feilchenfeld forceps



Code	Length
FML002	110 mm
FML003	90 mm



▶ “Kocher” forceps



Code	Length
FML004	140 mm

▶ Haemostatic forceps



- 1| They present a mechanism with teeth thanks to which once the forceps are closed they kept fixed
- 2| Once operation is finished they can be released easily
- 3| Haemostatic forceps are used to close blood vessels by applying pressure.

Code	Model	Length
FML005	Straight	130 mm
FML015	Curved	140 mm
FML017	Straight	140 mm

▶ Forceps for insects



Code	Length
FML006	100 mm

▶ Special forceps



CODE: FML007



CODE: FML010

Code	Description	Length
FML007	Spring effect	120 mm
FML010	Plain and wide tips	120 mm

▶ “Rochester-Pean” forceps



Code	Length
FML008	160 mm
FML020	140 mm



► Forceps with very thin tip




Code	Length
FML009	110 mm

► Maier clamp



Code	Length
FML018	180 mm

► Needle holders

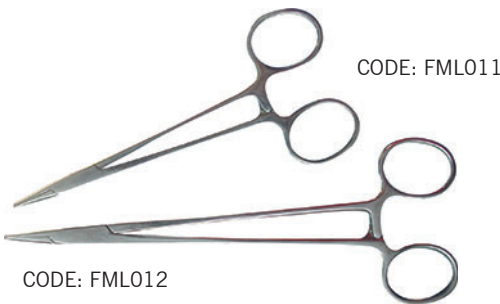

CODE: FML011

CODE: FML012

- 1 | They are used to take and hold surgical needles and are very similar to haemostatic forceps; the main difference is that needle holders present short and strong branches to hold the needle without damaging the suture material.
- 2 | The size of the needle holder must be in accordance with the needle size
- 3 | For this purpose, 2 models are available:
 - Baumgartner, for thin needles
 - Mayo-Hegar, for thick needles

Code	Length	Model
FML011	140 mm	Baumgartner
FML012	180 mm	Baumgartner
FML013	140 mm	Mayo-Hegar
FML014	180 mm	Mayo-Hegar
FML019	145 mm	Mayo-Hegar



CODE: FML013, FML014 and FML019