




Water distiller 4 l/h, QUARTZ

Code PDG002

- This is the water distiller with automatic and continuous production and quartz resistance that provides a high quality distillation. The results are **high purity water**, free of metal ions and pyrogenic substances able for preparation of reagents, bacteriological culture, final wash glassware, etc. The device is designed for clinical, medicine hospitals, researching laboratories, agro-alimentary, environmental or pharmaceutical industries and schools or universities.
- Is **made of high quality borosilicate glass** allowing continued use of the equipment and ensures its durability. It provides a **safety mechanism that automatically switches off the quartz resistance** when the cooling water pressure is reduced and reconnect when recover.
- Include an adapted bracket and small dimensions that allow **easy handling and cleaning**.

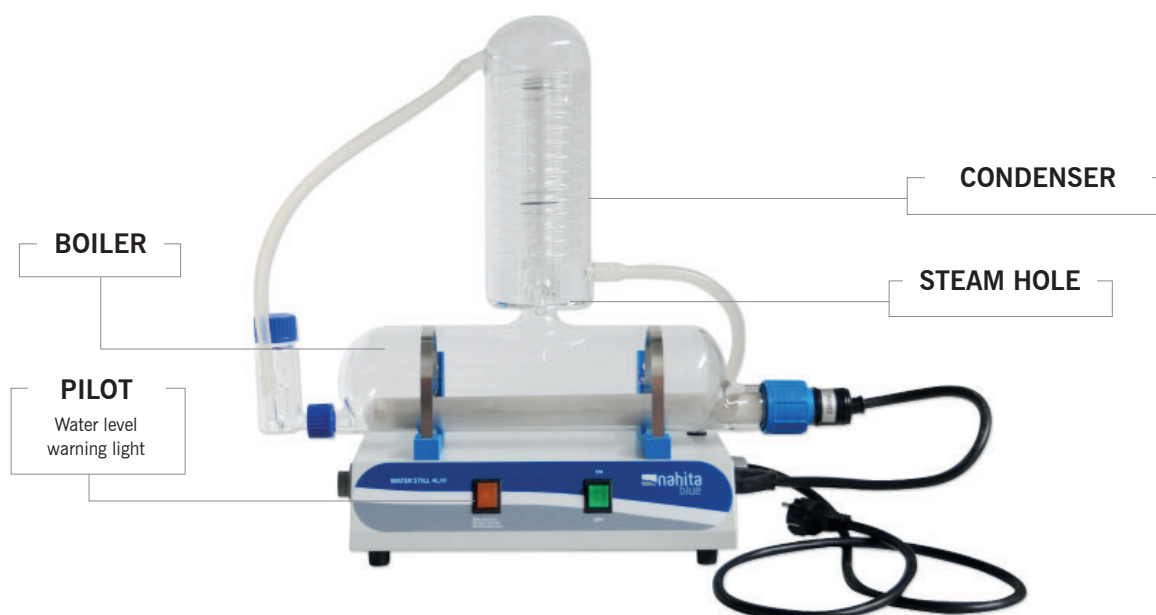

Technical specifications

Code	PDG002
Output	4 l/h
Conductivity at 20 °C	2 μ S/cm
Water consumption	1,0 litros/min
Resistance	1 x 2500 W
Power supply	230 V / 50-60Hz / 11 A
Dimensions	55 x 15 x 50 cm
Weight	4 Kg

Accessories

Referencia	Descripción
PDJ008	Boiler
PDJ009	Condenser
PDJ010	Level
PDJ011	Heater
PDJ012	Flow limiter
PDJ028	Guide bush for resistance for QUARTZ distiller
PDJ029	Plastic nut N. 42 for resistance, blue
PDJ030	O-ring for resistance, red for Quartz distiller