




Autoclave sterilization tape

Code ZCM008

The sterilisation indicator tape is specifically designed for use in autoclaves to monitor the sterilisation process in medical and laboratory environments.

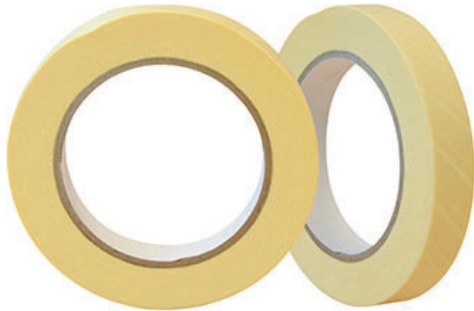

Characteristics :

- Tape manufactured with a heat-resistant adhesive that ensures adhesion to different surfaces.
- Colour change from white to brown when the appropriate sterilisation conditions are reached (autoclave 20 min at 121°C, 10 min at 125°C or 5 min at 135°).
- The indicator ink is lead-free, making it safe for both the environment and users.

Code	ZCM008
Total length (m)	50
Width (mm)	19