



## BUNSEN BURNERS

Code FFJ003, FFJ004, FFJ005, FFJ006, FFJ040

### Micro Bunsen Burner

Code FFJ003, FFJ004

- With adjustable air supply and a tapered serrated gas connection




Code	FFJ003	FFJ004
Gas type	Natural gas	GLP
Base material	Die-cast zinc alloy with metallic plating	
Inlet outer diameter	8 mm	
Thermal output power	800-1200 BTU	2000-3200 BTU
Weight	0.2 kg	
Height	89 mm	

### Adjustable Bunsen Burners

Code FFJ005, FFJ006

- With a threaded needle valve, aluminium mixing tubes for thorough air and gas mixing, and tapered serrated gas inlet connections for easy tube attachment

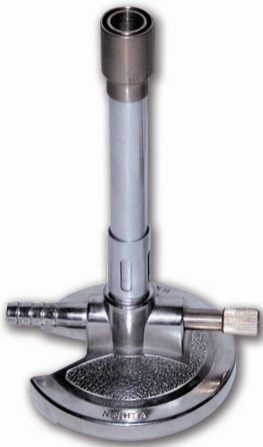

Code	FFJ005	FFJ006
Gas type	Natural gas	GLP
Base material	Die-cast zinc alloy with metallic plating	
Inlet outer diameter	8 mm	
Thermal output power	800-1200 BTU	2000-3200 BTU
Weight	0.3 kg	
Height	155 mm	

### Multigas Bunsen Burners with adjustable control

Code FFJ040

- With a threaded shut-off valve and flame stabilizer



Code	FFJ040
Gas type	All types of gas
Base material	Nickel-plated brass
Inlet outer diameter	1/2" & 1/4"
Thermal output power	300-3200 BTU
Weight	0.4 kg
Height	130 mm