



Bank of circular hotplates “Uniplac”

ADJUSTABLE TEMPERATURE HOTPLATE UP TO 400 °C.
SUITABLE FOR KJELDAHL DIGESTION ASSAYS AND
SOXHLET EXTRACTION.

FEATURES

Independent hotplate control using pulsed energy regulators.
Cast steel hotplates with embedded heating elements distributed throughout the hot-plate's surface. 
Adjustable height and angle stand bars for supporting condensers and other glass-ware etc.
External case made of AISI 304 stainless steel.

CONTROL PANEL

Independent temperature control for each position.
- Temperature control by pulsed energy regulators.
- Indicator lamps showing when the hotplate is on.

MODELS

Part No.	No. of places	Plate Ø cm	Height / Width / Depth (exterior) cm			Power W	Weight Kg
6002314	4	12	23	77	27	2200	15
6002316	6	12	23	105	27	3300	22



ACCESSORIES

Bath of AISI 304 stainless steel of 140 mm Ø x 75mm high, with concentric reduction ring lid, Suitable for water, oil or sand.
Part No. **7002510**



Round bottom flask adapter made of aluminium and suitable for Kjeldahl flasks of 50 to 100 ml and 50 to 100 ml round bottom flasks.
Part No. **1002290**



Round bottom flask adapter made of aluminium and suitable for Kjeldahl flasks of 300 to 800 ml and round bottom flasks of 250 to 500 ml.
Part No. **1002280**



ACCESSORIES FOR KJELDAHL DIGESTION

Kjeldahl support Bar for Kjeldahl flasks, made of AISI 304 stainless steel.
Part No. **6023141** 4 places
Part No. **6023161** 6 places



Clamp for Kjeldahl condenser tubes made of AISI 304 stainless steel.
Part No. **6000008**



Borosilicate glass fume tube.
Part No. **1000104** 4 place
Part No. **1000106** 6 place



EXTRACTION AND NEUTRALIZING SYSTEM FOR GASSES AND VAPOURS. Prevents water wastage.



“Scrubber”

Recirculating pump.

Features (see page 251).

Models

Unit “Scrubber”


Part No.	Height / Width / Depth (exterior) cm			Weight Kg
4001611	32	31	16	2

It is supplied with 3 Kg. acid vapours neutralizer solution.

Water recirculating vacuum pump

Part No.	Height / Width / Depth (exterior) cm			Vacuum level bar	Vacuum pump rate litres/minute	Weight Kg
4001612	44	39	28	0.98	10	10

ACCESSORIES SOXHLET EXTRACTION

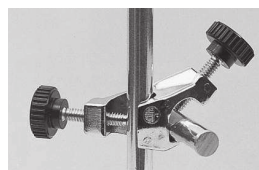

Boss head for supporting clamps.
Part No. **1000069**



Clamp for flasks or condenser tubes.
Part No. **7000403**