

COMECTA Motorized pipette dispenser "Rota-Filler 3000"®

FOR ANY KIND OF PIPETTES FROM 1 TO 100 ML.



FEATURES

Three positions switch for fast or slow aspiration and dispensation or gravitatorial dispensation for free dispensing.

The pressure over the ergonomic bulbs allows modulation, speed control of aspiration and dispensation.

It is supplied with a couple of 0,45 µm filter by default for an effective protection against contaminated air either of the equipment during aspiration or of the sample during dispensation.

Two rechargeable and interchangeable batteries.

Low-battery light indicator signal.

Can be used while recharging.



TECHNICAL DATA

Charger of 65/12 mA. Charge time: 14 hours.

Weight: 215 g.

Dimensions: 170 height x 47 width x 143 length mm.

Pipettes: from 1 to 100 ml (glass and plastic).

Motorized pipette dispenser "Rota-Filler 300"®.

Part no: **5900700**

COMPLEMENT

5 hydrophobe filters pack of 0,45 µm.

Part no: **5900701**



COMECTA Manual pipette dispenser "Pipet"

FOR ANY KIND OF PIPETTES FROM 1 TO 100 ML.

FEATURES

Ergonomic, comfortable, solid, lightweight and ease of use.

It has been designed with an ultracompressor silicone bulb that provides a soft manual control of the pipette, either in aspiration or dispensation.

Once the bulb is squeezed, control with the toggle both speeds, the aspiration and dispensation one. Push the middle button to expel blowing-out any liquid remaining in the pipette. 0,45 µm.

Measure: 25 mm Ø.


Weight: 130 g.

Manual pipette dispenser "Pipet". Part no **5900800**

COMPLEMENT

5 hydrophobe filters pack of 0,45 µm.

Part no: **5900801**

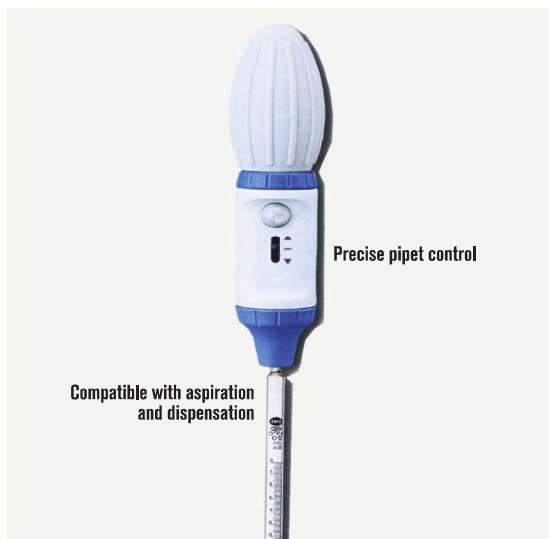

Measuring range "COM-PIP"

FEATURES

Made of polypropylene with silicon joints, for volumes up to 25 ml. Precise and safe, ergonomic design that allows filling and emptying without having to change the setting. Flexible thermoplastic seal against any type of burette. A button on the side empties the burette automatically. Easily dismantled and cleaned.

Part No.

Capacity

Identification colour

1001131

0 to 2 ml

Blue

1001132

0 to 10 ml

Green

1001133

0 to 25 ml

Red



PIPETTE RACK "COM-PIP"

Made from AISI 304 stainless steel.

Capacity 6 pipettes of differing sizes.

Part No. **1001702**

