




Graduated cylinders, PMP, Class A, tall form, red printed graduations

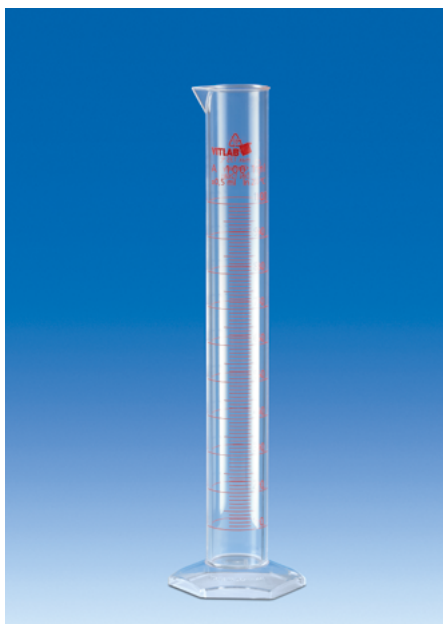

Highly transparent. DE-M marked.

With red printed graduations and ring marks at the primary scale points, calibrated 'In'.

The lot certificate supplied bears the batch number and the actual nominal value ascertained under the test conditions. The resulting deviations from the nominal value fall well under the allowed tolerances of Class A according to DIN 12681 and ISO 6706.

With printed batch number and year of production. Also available with DAKKS calibration certificate or individual quality certificate (at additional cost).

Hexagonal base with bottom studs provides high stability. To preserve markings, do not clean at temperatures exceeding 60 °C. Thus, conditionally autoclavable at 121 °C (2 bar) according to DIN EN 285. For autoclaving we recommend the design with molded graduations (Cat.-No. 64604 – 65304).

Volume ml	Tolerance ± ml	Graduation ml	Height mm	Ø mm	PU	Cat. No.
10	0.10	0.20	145	15	2	64614
25	0.25	0.50	170	22	2	64714
50	0.50	1.00	200	27	2	64814
100	0.50	1.00	250	33	2	64914
250	1.00	2.00	315	44	2	65014
500	2.50	5.00	360	58	2	65114
1000	5.00	10.00	440	69	1	65214
2000	10.00	20.00	535	97	1	65414



Graduated cylinders, PMP, Class A, tall form, molded graduations



Highly transparent. DE-M marked.

With molded graduations and ring marks at the primary scale points, calibrated 'In'.

The lot certificate supplied bears the batch number and the actual nominal value ascertained under the test conditions. The resulting deviations from the nominal value fall well under the allowed tolerances of Class A according to DIN 12681 and ISO 6706. With the laser engraved batch number and the year of manufacture. Also available with DAKKS calibration certificate or individual quality certificate (at additional cost).

Hexagonal base with bottom studs provides high stability. Thermal stress up to 121 °C (autoclaving) does not cause tolerance limits to be permanently exceeded.

Volume ml	Tolerance ± ml	Graduation ml	Height mm	Ø mm	PU	Cat. No.
10	0.10	0.20	145	15	2	64604
25	0.25	0.50	170	22	2	64704
50	0.50	1.00	200	27	2	64804
100	0.50	1.00	250	33	2	64904
250	1.00	2.00	315	44	2	65004
500	2.50	5.00	360	58	2	65104
1000	5.00	10.00	440	69	1	65204
2000	10.00	20.00	482	97	1	65304