


# DL-SUB01

## Horizontal Electrophoresis Tank

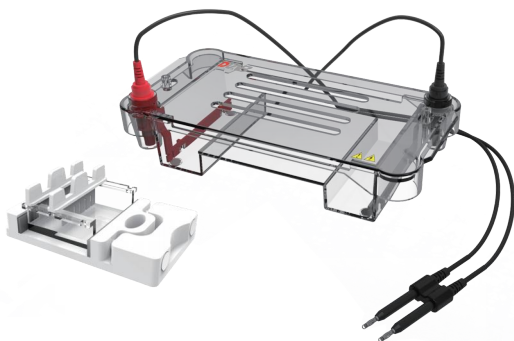

### Features

- Suitable for rapid separation, identification, preparation of DNA, and determination of its molecular weight
- Polycarbonate injection molding, no leakage
- The heat dissipation hole design on the upper cover effectively prevents condensation of water vapor
- Dedicated telescopic elastic glue maker, ensuring cleanliness and no glue leakage
- Limit function, effectively preventing reverse installation of positive and negative electrodes
- Detachable electrodes for easy replacement and maintenance
- 99.99% high-purity platinum electrode, the best conductivity, anti-aging

### Technical Parameters

Gel size: (W×L) : 60×75mm

Sample comb: 1.0/1.5mm 6/11wells

Buffer volume: 300mL

### Standard configuration

Name	Unit	Quantity
Lid with cables	set	1
Body tank with electrodes	set	1
Gel box	pc	1
Gel tray	pc	1
1.0mm thickness 6 wells comb	pc	2
1.0mm thickness 11 wells comb	pc	2
1.5mm thickness 6 wells comb	pc	2