



Glass wine pipette

Code From DLD001 to DLD008

Wine pipette is a glass instrument used in wineries to obtain samples of wine directly from the barrel for a later tasting or analysis. It consists of a long glass tube with a pair of handles at one end to hold the pipette while it is introduced into the barrel. At the opposite end of the pipette tip, the tube is closed except for a little hole over which the thumb is placed to avoid the pipette to empty when it is taken out the barrel.

Straight glass wine pipette

Code DLD001




Technical specifications

Code	Length	Capacity	Diameter
DLD001	640 mm	350 mL	30 mm

Straight glass wine pipette with rubber handle

Code DLD007, DLD008

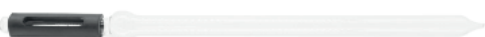

Technical specifications

Code	Length
DLD007	50 cm
DLD008	65 cm

Straight glass wine pipette with rubber handle, curved

Code DLD002, DLD003, DLD004, DLD005



Technical specifications

Code	Length
DLD002	25 cm
DLD003	35 cm
DLD004	45 cm
DLD005	65 cm